



**1** Choose proper size dressing after debridement.

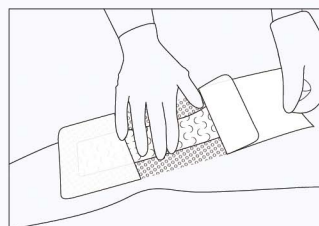

**2** Apply dressing directly onto the wound.

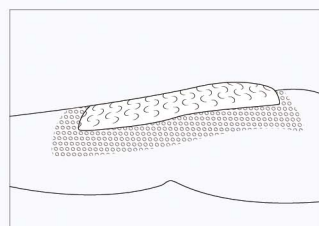

**3** Secure dressing with suitable bandage or tape if necessary.



Type	REF Code	Size (cm)	Size (in)
Adhesive	OP080060A	8 x 6 (5 x 3)	3.2" x 2.4" (2" x 1.2")
	OP100090A	10 x 9 (6 x 4.5)	4" x 3.6" (2.4" x 1.8")
	OP150090A	15 x 9 (10 x 4.5)	6" x 3.6" (4" x 1.8")
	OP100100A	10 x 10 (6 x 6)	4" x 4" (2.4" x 2.4")
	OP150100A	15 x 10 (10 x 6)	6" x 4" (4" x 2.4")
	OP200100A	20 x 10 (14.5 x 4.5)	8" x 4" (5.8" x 1.8")
	OP250100A	25 x 10 (19 x 4.5)	10" x 8" (7.6" x 1.8")
	OP300100A	30 x 10 (24 x 4.5)	12" x 8" (9.6" x 1.8")
	OP350100A	35 x 10 (29 x 4.5)	14" x 8" (11.6" x 1.8")
	OP175150A	17.5 x 15 (12.5 x 10)	7" x 6" (5" x 8")
	OP250200A	25 x 20 (19 x 14.5)	10" x 8" (7.6" x 5.8")
	OP500200A	50 x 20 (43 x 14.5)	20" x 8" (17.2" x 5.8")

Remark: Special specifications according to the order of the contract.

**FORYOU MEDICAL**  
**HUIZHOU FORYOU MEDICAL DEVICES CO., LTD**

HUIZHOU FORYOU MEDICAL DEVICES CO., LTD.  
ADD: No.1, Shangxia North Road, Dongjiang Hi-tech Industrial Park, Huizhou, Guangdong, China  
TEL: +86-752-5302185 FAX: +86-752-5302088 E-MAIL: medical2@foryougroup.com  
en.foryoumedical.com


**FORYOU MEDICAL**  
for you, for us, for life

**LUOFUCON®**

# Silicone Post-Op Dressing


**FORYOU MEDICAL**  
for you, for us, for life



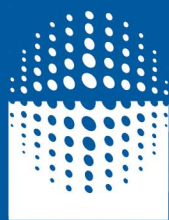

LUOFU CON®

# Silicone Post-Op Dressing

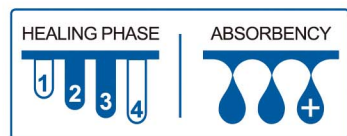

- Red Wound
- Yellow Wound
- Pink Wound
- M: Moisture Balance

## Flexible comfort for wound care after surgical operation.

LUOFU CON® Post-Operative Dressing is a soft, highly flexible and absorbent dressing. It is widely used to absorb and lock the blood and surgical exudates from surgical wounds. It's composed of a gentle silicone wound contact layer, non-woven wicking layer, a highly absorbent pad and a vapor permeable and water resistance film. The silicone wound contact layer is hypoallergenic and doesn't stick to the wound. It's very easy in application and cause no pain for dressing changes.

## Flexible

No move worries about the wounds on joints or difficult areas. High flexibility provides great comfort and adaptability.



## Crescent Texture

With high adaptability, the special design of crescent perforated pattern in the super absorbent pad creates dynamic flexibility as well as ensuring the dressing integrity during use. Horizontal and vertical absorption are also beneficial from this open-up structure.



### Protective Backing

High MVTR  
Smooth and flexible  
Water proof & bacteria barrier

### Super Absorbent Core

Exceptional absorption capability  
Helps minimize maceration risk  
Special crescent patterns provide dynamic flexibility

### Nonwoven Layer

Spread fluids to a wide surface  
Enhance absorption performance

### OptiSil Silicone Wound Contact Layer

Self-developed silicone gel and perforation technology  
Distributes wound exudate fast and evenly  
Gentle and secure with pain free removal  
Without Stripping of epidermal cells

## OptiSil

Self-developed silicone gel and perforation technology not only provides gentle care on the wound skin, but also enables exudate to run into the dressing amazingly quickly. No residues left after removal.

### Indications

LUOFU CON® Silicone Post-operative Dressing is indicated for management of moderate to high exuding wounds, including: leg ulcers, pressure ulcers, diabetic foot ulcers and dehisced surgical wounds.

### Benefits

- Excellent exudate management for Post-Operative wounds;
- Highly flexible and comfortable for joints;
- No second trauma for dressing changes;
- No sticking to the wounds or sutures;
- Easy for doctors to do wound inspection with one simple move.

### Contraindications

Do not use this product on dry or low exuding wounds, eyes, mucous membranes or in wound cavities because the dressing swells during absorption. Do not use on third-degree burns.

Should not be used on persons with a known allergy or hypersensitivity to the dressing or its components.